

How to stream sounds from your laptop to hearing aids using TV streamer (RCA to Stereo converter)



Step 1

Step 2


Step 3

Step 4

Step 5

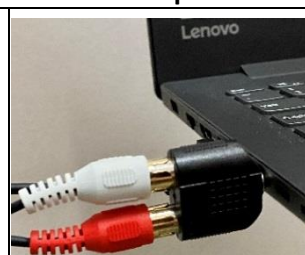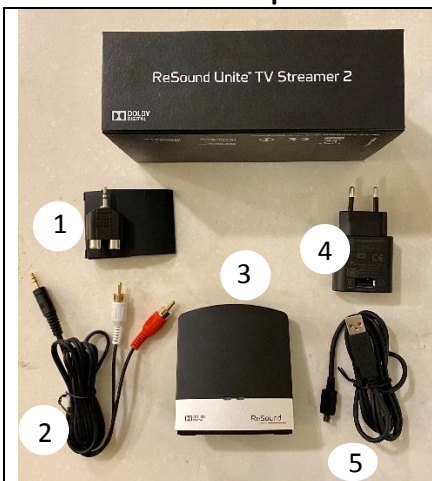

Close up

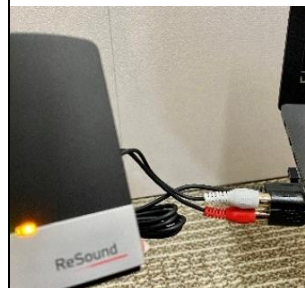

What you find in Resound TV streamer box for the setup:

- 1) RCA-to-Stereo converter
- 2) RCA cable (red and white)
- 3) TV streamer
- 4) Power plug
- 5) Cable for power plug

Connect the stereo port end of the RCA cable to the TV streamer's audio port located at the back of the TV streamer (as indicated by the arrow)

Connect the RCA's red and white ports to the RCA-to-Stereo converter

Connect the stereo port of the RCA-to-Stereo converter to your laptop's audio port (marked under (6))

This audio port is also known as the 3.5mm headphone jack

You need to power up the TV streamer next.

You can either connect to the laptop using the usb cable marked as (5) or use the 2-pin power plug marked under (4) outlined in step 1 to connect to the wall's pocket socket

The TV streamer will light up and ready to stream sounds to your hearing aids once this is done successfully. Enjoy the new listening experience!